


ERGO Ceramic

Ceramic modeling instrument

Double-sided special instruments for modelling porcelain with three replaceable instrument tips.

Advantages

- Instrument No. 1: Large blade for separating and reducing (replaceable). A large claw scraper for precise resection when modelling the tooth shape.
- Instrument No. 2: Fine nerve needle for creating fissures and enamel cracks (replaceable). Small blade for working on areas which are difficult to reach (replaceable).

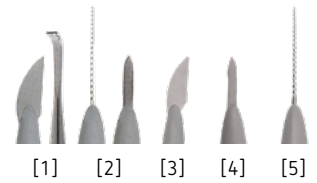

Ordering information

| | | |
|------------------------------------|--------------------|--------------|
| ERGO Ceramic set | with 2 instruments | No. 11611000 |
| [1] ERGO Ceramic Instrument No. 1 | 1 piece | No. 11611100 |
| [2] ERGO Ceramic Instrument No. 2 | 1 piece | No. 11611200 |
| [3] Spring tip „large blade“ | 2 pieces | No. 11611110 |
| [4] Spring tip „small blade“ | 2 pieces | No. 11611210 |
| [5] Spring tip „fine nerve needle“ | 2 pieces | No. 11611220 |

ERGO Acryl

Instrument for veneering composite

Double-sided modelling instruments ideal for acrylic facings.

Advantages

- Three individual flexible blades, 0.2 mm thick.
- Flexible blades are secured to the instrument handle by a threaded connection and can be easily exchanged, as required.
- Ideal for applying and modelling modern facing acrylics.



Ordering information

| | | |
|----------------------------------|----------|--------------|
| ERGO Acryl set | | No. 10521000 |
| [1] ERGO Acryl Instrument, No. 1 | 1 piece | No. 10521100 |
| [2] ERGO Acryl Instrument, No. 2 | 1 piece | No. 10521200 |
| [3] ERGO Acryl Instrument, No. 3 | 1 piece | No. 10521300 |
| [4] Flexible blade tip, No. 1 | 2 pieces | No. 10521110 |
| [5] Flexible blade tip, No. 2 | 2 pieces | No. 10521210 |
| [6] Flexible blade tip, No. 3 | 2 pieces | No. 10521310 |


ERGO Wax

Wax modeling instrument

ERGO wax instruments are characterised by the wide variety of applications for which they can be employed and by their individual adaptability.

Advantages

- Instruments can be ground according to individual needs.
- Probes, blades and scrapers made from handground stainless steel.
- Fine, anodised handles, comfortable to hold, heat-resistant and insulated.

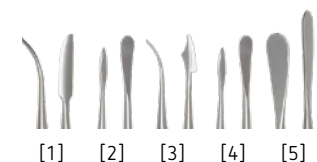

Ordering information

| | | |
|-------------------------------|---------|--------------|
| [1] ERGO Wax Instrument No. 1 | 1 piece | No. 10342001 |
| [2] ERGO Wax Instrument No. 2 | 1 piece | No. 10342002 |
| [3] ERGO Wax Instrument No. 3 | 1 piece | No. 10342003 |
| [4] ERGO Wax Instrument No. 4 | 1 piece | No. 10342004 |
| [5] ERGO Wax Instrument No. 5 | 1 piece | No. 10342005 |